



ApoE genotyping in DNA samples from PPMI

Dena G Hernandez

National Institute on Aging, Laboratory of Neurogenetics

Summary

ApoE genotyping was performed on 304 DNA samples collected according to the PPMI Research Biomarkers Laboratory Manual. Genotypes were determined with the use of allele-specific oligonucleotide probes labeled with fluorogenic reporter dyes (TaqMan method).

Method

Two non-synonymous single nucleotide polymorphisms (SNPs), rs429358 (*APOE*-C112R) and rs7412 (*APOE*-R158C) were genotyped in each sample in order to distinguish between *APOE* ϵ 2, ϵ 3 and ϵ 4 alleles. Taqman Assays (Applied Biosystems Assay-On-Demand part numbers C__3084793_20 and C__904973_10) were used per manufacturers protocol to genotype these SNPs on a 7900HT Sequence Detection System (Applied Biosystems). PCR amplification, plate read and allelic discrimination were performed on SDS instrumentation using the Computer Software SDS V2.4 2010.


References

Livak KJ. Allelic discrimination using fluorogenic probes and the 5' nuclease assay. *Genet Anal.* 1999 Feb;14(5-6):143-9. PubMed PMID: 10084106.

Taqman SNP Genotyping Assays Protocol:

<http://www3.appliedbiosystems.com/sup/URLRedirect/index.htm?xDoD=4332856>

Parkinson's Progression Marker Initiative Research Biomarkers Laboratory Manual (Biologic Manual) <http://www.ppmi-info.org/wp-content/uploads/2011/05/PPMI-Biologics-Manual-April-2011-FINAL.pdf>

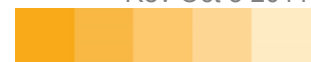



PARKINSON'S
PROGRESSION
MARKERS
INITIATIVE

Play a Part in Parkinson's Research

About the Authors

This document was prepared by Dena G Hernandez, National Institute on Aging, Laboratory of Neurogenetics. For more information please contact Dena Hernandez by email at

hernand@mail.nih.gov

Notice: This document is presented by the author(s) as a service to PPMI data users. However, users should be aware that no formal review process has vetted this document and that PPMI cannot guarantee the accuracy or

